

# BODYPOINT®



# PRODUCT GUIDE

# BODYPOINT® PRODUCT GUIDE



## 4-9

### PELVIC POSITIONING

- Hip Belts
- Sub-ASIS Pads and Belts
- Slip-On Pads
- Hip Belt Accessories
- Pivot Mount Belts
- Hip Grip™
- Leg Harnesses

## 10-13

### LOWER BODY POSITIONING

- Ankle Huggers®
- Adjustable Angle Footplates
- Calf Panels and Calf Straps

## 14-18

### UPPER BODY POSITIONING

- Stayflex™ Trunk Supports
- Shoulder Harnesses
- Chest Belts
- Monoflex™
- Aeromesh™ Shower Chair Belts

## 19-21

### POWERCHAIR COMPONENTS

- Midline Joystick Mounting System
- Joystick Handles


## 22

### HARDWARE

- Fixits™ and Clickits®
- FrameSaver™ Clamps
- Frame Mounting Clamps
- Quick-Release Connectors

## 23

### TERMS AND CONDITIONS

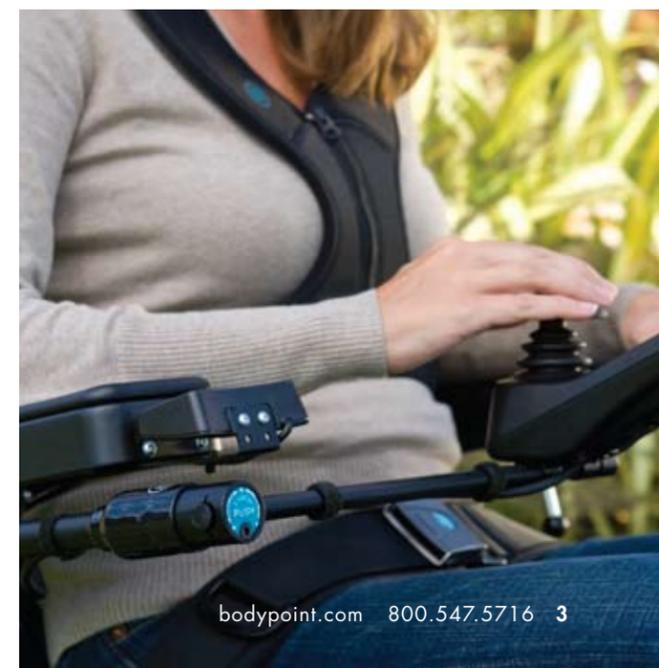

## PELVIC POSITIONING

The pelvis is the key to postural alignment as it dictates the position of the trunk, head and extremities.

## PUSH-BUTTON BUCKLE

The new Bodypoint® Push-Button Buckle is sleek and low profile. Engineered for comfort and performance, this new design offers a range of features that therapists and users have requested for years. The patented construction allows for removal of the buckle cover for easy cleaning, or to change the desired button hole size—in effect providing three buckle options in a single design.

Light Weight • Durable • Low Profile • Easy to Clean • Assisted Spring Release



Lighter weight (1/3 lighter than comparable push-button models) with rounded edges that improve comfort and prevent snagging on clothing.



Chrome-plated steel construction for superior strength and durability. Lower profile for a less bulky appearance.



Optional covers: Reduced Access Cover (Ø 14mm button), Security Cover (Ø 6mm button). Sold separately. See page 9 for ordering information.

# SELECT A BODYPOINT HIP BELT IN 6 STEPS

## 1. SELECT BELT TYPE

- Padded (pg. 7)
- Non-Padded (pg. 8)
- Hook-and-Loop-Compatible (pg. 8)

## 2. TWO OR FOUR-POINT



### TWO-POINT HIP BELT

Two points of attachment to the wheelchair. Suitable when the belt is supporting users with low tone or weakness.



### FOUR-POINT HIP BELT

Four points of attachment to the wheelchair. Primary straps position the belt like two-point hip belt and secondary straps anchor the belt in position to prevent it from riding up into the abdomen or twisting.

## 3. SELECT BUCKLE

OPTIONAL COVERS:



### Push-Button Buckle

Low profile, smooth, snag-free corners. Removable cover for changing release hole size and easy cleaning. Standard Access Cover (a) included. Reduced Access Cover (b) and Security Cover (c) sold separately.



### Plastic Side-Release

Lightweight, secure snap action.



### Rehab Latch™ Buckle

Lightweight, strong, low profile, easy to operate with limited hand function.



### Extra Small Metal Push-Button

Strong, secure, child-sized buckle. Buckle Security Cover also available. (See pg.9)



### Compact Push-Button

Firm button action prevents accidental release.

## 4. SELECT A PULL TYPE



### CENTER-PULL HIP BELT

**Operator:** User or caregiver

**Adjustment:** Tightens at the center buckle.

**Buckles:** Available in the full buckle range.

**Applications:** Clients with low muscle tone or weakness.



### REAR-PULL HIP BELT

**Operator:** Caregiver

**Adjustment:** Tightens from the rear of the pad pulling toward the center.

**Buckles:** Available in full buckle range except for Extra Small Metal Push-Button.

**Applications:** Clients with high tone. Counters pelvic rotation by allowing the caregiver to position one side of the pelvis, lock it in place and then position on the other side.



### DUAL-PULL HIP BELT

**Operator:** User or caregiver

**Adjustment:** Tightens from both sides of the center buckle.

**Buckles:** Side-release buckle only

**Applications:** Clients requiring greater adjustment to accommodate clothing and weight changes.

## 5. SELECT A SIZE

Sizes available from extra small to large. See measurement section to the left for details.

## 6. SELECT A MOUNTING OPTION



**Cinch-Mount™**  
Time-saving, precise, easy adjustment

**Flat-Mount™**  
Versatile, lowest profile



### FrameSaver™ Clamp (optional)

Create a mounting point without drilling. Fits Ø7/8"-Ø1" (22-25mm) tube.

## TWO-POINT PADDED HIP BELTS

		Push-Button	Plastic Side-Release	Rehab Latch™	Extra Small Metal Push-Button	Compact Push-Button
CENTER-PULL	XS		HB215XS		HB235XS <sup>1</sup>	
	S	HB205S*	HB215S	HB225S		HB245S
	M	HB205M*	HB215M	HB225M		HB245M
	L	HB205L**	HB215L	HB225L		HB245L
DUAL-PULL	XS		HB216XS			
	S		HB216S			
	M		HB216M			
	L		HB216L			
REAR-PULL	XS		HB219XS		HB239XS <sup>1</sup>	
	S	HB209S*	HB219S	HB229S		HB249S
	M	HB209M*	HB219M	HB229M		HB249M

\*Accepts the FS018D-1 Reduced Access and FS018C-1 Security Covers (see page 9). \*\*Accepts the FS019D-1 Reduced Access and FS019C-1 Security Covers (see page 9). <sup>1</sup>Accepts FS029C-1 Buckle Security Cover (see page 9).

## FOUR-POINT PADDED HIP BELTS

		Push-Button	Plastic Side-Release	Rehab Latch™	Extra Small Metal Push-Button	Compact Push-Button
CENTER-PULL	XS		HB415XS		HB435XS <sup>1</sup>	
	S	HB405S*	HB415S	HB425S		HB445S
	M	HB405M*	HB415M	HB425M		HB445M
	L	HB405L**	HB415L	HB425L		HB445L
DUAL-PULL	XS		HB416XS			
	S		HB416S			
	M		HB416M			
	L		HB416L			
REAR-PULL	XS		HB419XS		HB439XS <sup>1</sup>	
	S	HB409S*	HB419S	HB429S		HB449S
	M	HB409M*	HB419M	HB429M		HB449M

\*Accepts the FS018D-1 Reduced Access and FS018C-1 Security Covers (see page 9). \*\*Accepts the FS019D-1 Reduced Access and FS019C-1 Security Covers (see page 9). <sup>1</sup>Accepts FS029C-1 Buckle Security Cover (see page 9).

# PADDED HIP BELTS


To complete your part number, select the end-fitting option below:

- Flat-Mount + FrameSaver Clamp + -A1
- Flat-Mount Only -B1
- Cinch-Mount + FrameSaver Clamp + -A2
- Cinch-Mount Only -B2

Example: **HB205S-B1**

Product number for the Center-Pull Push-Button Buckle, size Small with Flat-Mount only.

# HOOK-AND-LOOP-COMPATIBLE BELTS

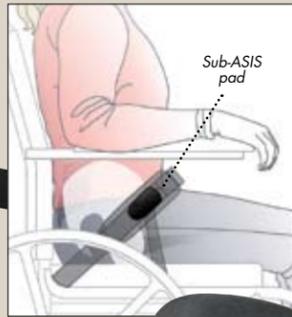

**SUB-ASIS PADS**  
The teardrop-shaped gel pads comfortably space the belt away from the abdomen (reducing pressure on the bladder), while firmly controlling the pelvis just below the ASIS.

PD402-2 Extra Small/Small  
PD401-2 Medium/Large

Combined with Sub-ASIS pads, our hook-and-loop-compatible belt provides two firm points of contact with the pelvis. The result is greatly enhanced stability and improved control of obliquity and rotation while minimizing abdominal compression.

		Size	Push-Button	Plastic Side-Release	Rehab Latch™	Extra Small Metal Push-Button
TWO-POINT	Center-Pull	M	LB205M*		LB225M	
		L			LB225L	
	Rear-Pull	M			LB229M	
FOUR-POINT	Center-Pull	XS				LB435XS¹
		M	LB405M*		LB425M	
	L	LB405L**		LB425L		
	Dual-Pull	XS		LB416XS		
	Rear-Pull	M			LB429M	

\*Accepts the FS018D-1 Reduced Access and FS018C-1 Security Covers (see page 9). \*\*Accepts the FS019D-1 Reduced Access and FS019C-1 Security Covers (see page 9). ¹Accepts FS029C-1 Buckle Security Cover (see page 9).

To complete your part number, select the end-fitting option below:

Flat-Mount + FrameSaver Clamp		-A1
Flat-Mount Only		-B1
Cinch-Mount + FrameSaver Clamp		-A2
Cinch-Mount Only		-B2

# NON-PADDED HIP BELTS

To complete your part number, select the end-fitting option in the chart above



**SLIP-ON PADS**  
Allow you to add cushioning where additional relief is needed on belts and harnesses.

PD304 Small  
Fits webbing: 1"(25mm)  
Length: 8"(20cm)  
PD305 Medium  
Fits webbing: 1-1/2"(38mm)  
Length: 10"(25cm)  
PD306 Large  
Fits webbing: 2"(50mm)  
Length: 12"(30cm)

	Size	Push-Button	Plastic Side-Release	Rehab Latch™	Compact Push-Button
Center-Pull	XS		NB215XS		
	M	NB205M*	NB215M	NB225M	NB245M
	L	NB205L**	NB215L	NB225L	NB245L¹
Center-Pull Quad Belt	L			NB225L-T	

\*Accepts the FS018D-1 Reduced Access and FS018C-1 Security Covers (see page 9). \*\*Accepts the FS019D-1 Reduced Access and FS019C-1 Security Covers. (See page 9)

NEW

**OPTIONAL PUSH-BUTTON BUCKLE COVERS**  
(Push-Button Buckle shown on page 6)



**REDUCED ACCESS COVER**  
FS018D-1 1-1/2", Ø14mm button  
For use with small and medium belts.  
FS018D-50 Pack of 50  
FS019D-1 2", Ø14mm button  
For use with large belts.  
FS019D-50 Pack of 50



**SECURITY COVER**  
FS018C-1 1-1/2", Ø6mm button  
For use with small and medium belts.  
FS018C-50 Pack of 50  
FS019C-1 2", Ø6mm button  
For use with large belts.  
FS019C-50 Pack of 50



**EXTRA SMALL BUCKLE SECURITY COVER**  
Snaps securely onto our Extra Small Metal Push-Button Buckle, limiting access by even the smallest fingers. Open instantly with a pen tip. FS029C-1



Shown: Previous Metal Push-Button Buckle with Buckle Security Cover Attached.



**BUCKLE SECURITY COVER**  
Snaps securely onto our previous Metal Push-Button Buckle (no longer available), limiting access by even the smallest fingers. Open instantly with a pen tip. FS025C-1

# PIVOT-MOUNT BELTS

Provides direct mounting to rigid surfaces, requiring only 1/4" (6mm) hole.



Rehab Latch™  
SB225M Medium  
SB225L Large  
SB225XL Extra Large




Extra Small Metal Push-Button  
SB235XS Extra Small

# REACH BEYOND THE CHAIR with the Bodypoint® Hip Grip™



US and International patents and patents pending

**HIP GRIP**  
An innovative dynamic pelvic stabilization device to help the wheelchair user maintain pelvic stability while allowing functional pelvic movements. The Hip Grip allows the pelvis to perform a natural anterior and posterior range of motion, while providing variable resistance to return the pelvis to its neutral posture. The combination of improved pelvic positioning and dynamic pelvic movement enhances sitting posture and improves functional activities, such as reaching and propelling the wheelchair. For additional information, contact Bodypoint.

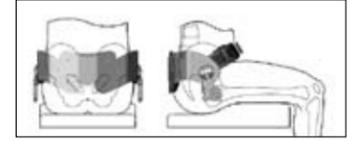

# LEG HARNESS

Holds each leg at the inner thigh, providing exceptional resistance to pelvic rotation and obliquity, without limiting anterior movement or pressing on the bladder. Sold in pairs.

LH600S Small 9" (23cm)  
LH600M Medium 13" (33cm)  
LH600L Large 15" (38cm)



# MEASURING FOR A LEG HARNESS




Measure from the Anterior Superior Iliac Spine (ASIS) to the point where the inner thigh touches the seat.

S Fits 8-11" (20-28cm)  
M Fits 11-14" (28-36cm)  
L Fits 14-16" (36-41cm)

## LOWER BODY POSITIONING

Lower body positioning affects the position of the pelvis, balance, and freedom of movement of the upper body.

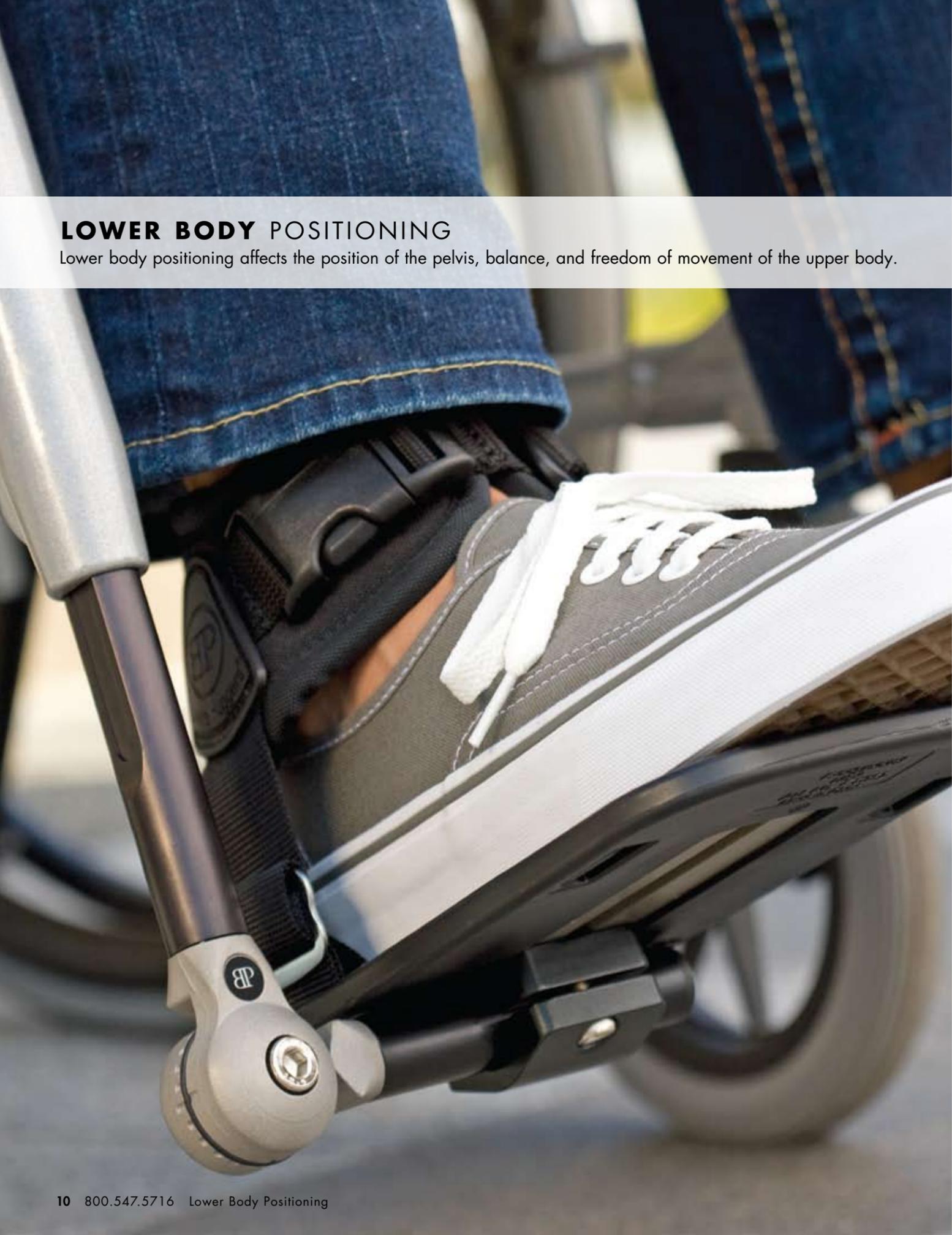

## ANKLE HUGGERS®

Ankle Huggers are the original patented solution for stabilizing the feet while allowing controlled movement. They reduce joint stress by holding above the ankle, and greatly improve comfort and tolerance for positioning. Engineered for outstanding fit and lasting quality with our patented Laminar Pad™ and proprietary webbing. Ankle Huggers are the most secure means available for holding the foot.



Patented reinforcement patch provides enhanced strength between the straps and pads.



Circumferential pad shape provides even pressure distribution.



Side-release attachment for exact fit. Also available in hook-and-loop closure for easy adjustment.



### HOOK-AND-LOOP RELEASE ANKLE HUGGERS

Sold in pairs (1 right/1 left)  
Soft hook-and-loop adjustment strap conforms to ankle curvature and is quickly tightened with each use to ensure a snug fit.

- FT210XS *Extra Small*
- FT210S *Small*
- FT210M *Medium*
- FT210L *Large*
- FT210XL *Extra Large*



### SIDE-RELEASE ANKLE HUGGERS

Sold in pairs (1 right/1 left)  
Adjustment strap retains its precise position with every use, and side-release buckle ensures strong closure even with extensive use.

- FT220XS *Extra Small*
- FT220S *Small*
- FT220M *Medium*
- FT220L *Large*
- FT220XL *Extra Large*



### AEROMESH™ CALF PANEL

Features our proprietary Aeromesh material. Aeromesh fabric is bonded to a durable stiffener to allow air flow while maintaining shape.

- SP103S *Small*
- SP103M *Medium*
- SP103L *Large*



### AEROMESH PADDED CALF STRAPS

Calf Straps provide comfortable support, keeping feet safely on the footplates. Retains its shape for the best pressure distribution and durability.

- SP102S *Small*
- SP102M *Medium*
- SP102L *Large*



### HOOK-AND-LOOP TOE STRAPS

Toe strap slots are provided on all Bodypoint footplates to allow for easy attachment. Sold in pairs.

FT221

## MEASURING FOR ANKLE HUGGERS®



Measure circumference above ankle.

XS	5-1/2"-7" (14cm-17cm)
S	6-1/2"-8" (17cm-20cm)
M	7-1/2"-9" (19cm-23cm)
L	8-1/2"-11-1/2" (22cm-29cm)
XL	11"-13" (28cm-33cm)

## FOOTPLATE MOUNTING OPTIONS



Mounted in hanger tubes



Mounted high for short leg lengths (requires HW404-2 or HW406-2 clamp kits)



Mounted at different heights for single leg amputee (requires HW404-2 or HW406-2 clamp kits)



### L-SHAPED EXTENSION TUBE

Ø3/4" (19mm)

- FT021 Length: 8-1/2" (22cm)
- FT022 Length: 12-1/2" (32cm)



### QUICKIE EXTENSION TUBE

Ø3/4" (19mm)

- FT026 Length: 12-1/2" (32cm)



### FULCRUM SERIES EXTENSION TUBE

Ø3/4" (19mm)

- FT023 Length: 8-1/2" (22cm)
- FT024 Length: 12-1/2" (32cm)



### MOUNTING CLAMP

HW404 Ø3/4" (19mm) X Ø7/8" (22mm)  
HW404-2 Pair

HW406 Ø3/4" (19mm) X 1" (25mm)  
HW406-2 Pair



### FOOTPLATE CLAMP

Clamp permits easy adjustment of footplate and locks tightly in position for absolute stability.

FT032 Ø3/4" (19mm)



### FOOTPLATE

Heavy gauge aluminum with super-tough, non-skid nylon surface stands up to hard use.

- FT051 *Small*  
Width: 4" (10cm)  
Length: 7-3/4" (20cm)
- FT052 *Large*  
Width: 4-3/4" (12cm)  
Length: 9-1/4" (24cm)

## ADJUSTABLE-ANGLE FOOTPLATE KITS

Includes: choice of extension tubes (pair), choice of footplates (pair), and footplate clamps (pair).


L-SHAPED		FULCRUM, Ø3/4" (19MM) TUBE		FULCRUM, QUICKIE	
Footplate Size	Tube Length	Footplate Size	Tube Length	Footplate Size	Tube Length
	8-1/2" (22cm) 12-1/2" (32cm)		8-1/2" (22cm) 12-1/2" (32cm)		12-1/2" (32cm)
S	FT101 FT101L	S	FT105 FT105L	S	FT105LQ
L	FT102 FT102L	L	FT106 FT106L	L	FT106LQ



### ADJUSTABLE ANGLE FOOTPLATE


- FT111 *Small*
- FT112 *Large*

## MEASURING FOR CALF PANELS AND CALF STRAPS

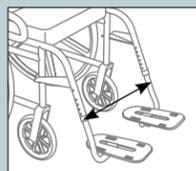


For proper fit, measure chair width.

S	12"-15" (31cm-38cm)
M	15"-19" (38cm-48cm)
L	19"-22" (48cm-56cm)

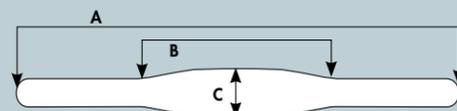

Calf Strap Dimensions:

	A	B	C
S	31-1/2" (80cm)	14" (36mm)	3-1/2" (89mm)
M	35-1/2" (90cm)	16" (41mm)	3-1/2" (89mm)
L	39-1/2" (100cm)	18" (46mm)	3-1/2" (89mm)

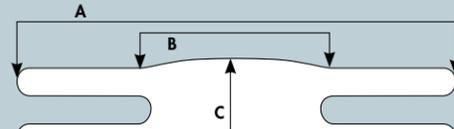

Calf Panel Dimensions:

	A	B	C
S	31-1/2" (80cm)	14" (36mm)	7-1/2" (19cm)
M	35-1/2" (90cm)	16" (41mm)	8-1/4" (21cm)
L	39-1/2" (100cm)	18" (46mm)	9-1/4" (23cm)



## UPPER BODY POSITIONING

Trunk support is essential for good respiration, head positioning, and functional reach with stability.



### STAYFLEX™ TRUNK SUPPORT

The patented design of the Bodypoint® Stayflex includes a stabilizing lower panel to limit stretching and riding up. The proprietary elastic pad material allows a comfortable, controlled amount of stretch. Plus, the lower quick-release swivel buckles self-adjust to the line of tension for proper fit and placement.



Proprietary elastic pad material offers controlled stretch for dynamic movement.



Rear-pull attachment is ideal for use with thicker backrests or when you want the buckle or adjustment away from the user.




Stabilizing lower panel limits stretching and riding up.

STANDARD STAYFLEX	FRONT-PULL		REAR-PULL		
	Size	With Zipper	Without Zipper	With Zipper	Without Zipper
 <i>Shown with zipper</i>	XS		SH311XS		SH321XS
	S	SH310S	SH311S	SH320S	SH321S
	M	SH310M	SH311M	SH320M	SH321M
	L	SH310L	SH311L	SH320L	SH321L
	XL	SH310XL	SH311XL	SH320XL	SH321XL

NARROW STAYFLEX	FRONT-PULL		REAR-PULL		
	Size	With Zipper	Without Zipper	With Zipper	Without Zipper
 <i>Shown with zipper</i>	M	SH330M	SH331M	SH340M	SH341M
	L	SH330L	SH331L	SH340L	SH341L
	XL	SH330XL	SH331XL	SH340XL	SH341XL

### MEASURING FOR STAYFLEX™ AND SHOULDER HARNESS

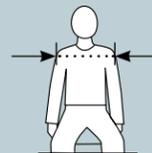

Measure shoulder width. If the shoulder width falls between two sizes, consider other factors such as growth, weight changes, and clothing.

XS	9-1/2-11" (24-28cm)
S	11-13" (28-33cm)
M	13-16" (33-41cm)
L	16-19" (41-48cm)
XL	19-22" (48-56cm)

### STAYFLEX™ AND SHOULDER HARNESS ATTACHMENT OPTIONS

NEW REAR-PULL OPTION



**REAR-PULL**  
Rear-pull attachment is ideal for use with thicker back-rests or when you want the buckle or adjustment away from the user.



**OPTIONAL CONFIGURATION**  
The included side-release buckles may be removed from the rear-pull strap to provide a continuous length of straight webbing for use with cam buckles or other mounting hardware.

Cam buckles shown not included. See page 22 for ordering information.



**FRONT-PULL**  
Front-pull attachment provides the maximum stability while allowing adjustment standing either behind or in front of the user.

STANDARD H-STYLE SHOULDER HARNESS	Size	Pad Length	Front-Pull	Rear-Pull
	S	12" (31cm)	SH210S	SH220S
M	13-1/2" (34cm)	SH210M	SH220M	
L	16" (41cm)	SH210L	SH220L	
XL	18-1/2" (47cm)	SH210XL	SH220XL	

Has an adjustable sternum strap to ensure a perfect fit; ideal for controlling strong, spastic movement.

**TRIMLINE STYLE SHOULDER HARNESS** Made of two independent padded straps which allow the greatest degree of freedom in fitting the chest.


TRIMLINE STYLE SHOULDER HARNESS	Size	Pad Length	Front-Pull	Rear-Pull
	S	18-1/2" (47cm)	SH230S	SH240S
M	11-1/2" (29cm)	SH230M	SH240M	
L	13-1/2" (34cm)	SH230L	SH240L	



**CROSSOVER POSITION**  
Comfortable fit which does not restrict arm movement.


**BACKPACK POSITION**  
Promotes respiration by encouraging shoulder retraction while stabilizing the trunk.

### STAYFLEX MEASUREMENTS

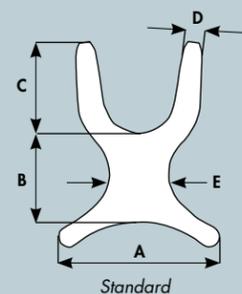

#### STANDARD STAYFLEX

Size	A	B	C	D	E
XS	1-1/4" (22cm)	5 (13mm)	5-1/4 (13mm)	1-1/2 (4mm)	3-3/4 (9mm)
S	9" (23cm)	5-1/2" (14cm)	6-1/2" (16cm)	1-3/4" (4cm)	4-3/4" (12cm)
M	12" (30cm)	7" (18cm)	8" (20cm)	1-3/4" (4cm)	5" (13cm)
L	15" (38cm)	8-1/2" (21cm)	9" (23cm)	2-1/4" (6cm)	5-1/4" (13cm)
XL	18" (46cm)	10" (25cm)	10-1/2" (27cm)	2-1/4" (6cm)	6-1/2" (16cm)

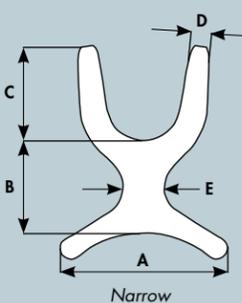

#### NARROW STAYFLEX

Size	A	B	C	D	E
M	12" (30.5cm)	7" (18cm)	8" (20.5cm)	1-3/4" (4cm)	3-1/2" (9cm)
L	15" (38cm)	8-1/2" (21.5cm)	9" (23cm)	2-1/4" (6cm)	3-1/2" (10cm)
XL	18" (46cm)	10" (25.5cm)	10-1/2" (27cm)	2-1/4" (6cm)	3-3/4" (10cm)

### FEEL SAFE AND SECURE with the Bodypoint® Shower Chair Belt




#### SHOWER CHAIR BELT

The Bodypoint shower chair belt with our proprietary Aeromesh™ padding provides safety without compromising the comfort and health of the user. Water and air easily pass through the Aeromesh fabric, allowing the client to be bathed completely and the belt to dry quickly after use. Its flexible construction and universal shape allow it to be used in various ways to secure the user: across the chest or abdomen to prevent forward/lateral falls, or fitted snugly across the knees to keep clients from slipping forward and out of the chair.



- BB101M-22-1 Medium (22mm clamps)
- BB101M-25-1 Medium (25mm clamps)
- BB101L-22-1 Large (22mm clamps)
- BB101L-25-1 Large (25mm clamps)
- BB101M-1 Medium (without clamps)
- BB101L-1 Large (without clamps)

Water and air easily pass through the Aeromesh™ fabric.

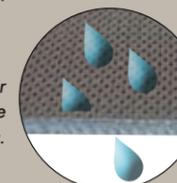



### MONOFLEX ANTERIOR TRUNK SUPPORT

The perfect dynamic abdominal or chest belt. Comfortable, high-strength elastic pad provides just the right amount of stretch while maintaining substantial control. Side-release buckles allow adjustment and release on both sides for convenience. Versatile mounting hardware fits any type of solid or sling backrest.

Size	Pad Length	Pad Width	Part No.
XS	12-1/4" (31cm)	3" (76mm)	SH103XS
S	15-1/2" (39cm)	3-1/2" (89mm)	SH103S
M	18-1/4" (46cm)	4" (10cm)	SH103M
L	21-1/2" (55cm)	4-1/2" (12cm)	SH103L
XL	24-1/2" (62cm)	5" (13cm)	SH103XL

### MEASURING FOR MONOFLEX™

Measure the user's trunk around the front and sides at the height at which the Monoflex will be worn. Consult the pad sizes (above) and select the size that will cover the front and wrap at least halfway around each side of the trunk, leaving a gap in the front of the wheelchair's backrest to allow space for tightening the buckles.



Size	Min/Max Fit Range
XS	114- 17-1/2" (36-45cm)
S	18- 21-1/2" (46- 55cm)
M	22- 25-1/2" (56-65cm)
L	26- 29-1/2" (66-75cm)
XL	30- 33-1/2" (76-85cm)



### ELASTIC CHEST STRAP

A versatile, easily removable anterior trunk support, abdominal support or positioning tool for therapists, users and caregivers. The simple thumb loop hook-and-loop closure provides easy one-handed opening and adjustment. Usable Length: 56" (152cm)

- SP110S** Small Width: 3" (75mm)
- SP110L** Large Width: 4-1/2" (11cm)



### CHEST BELT

Features a premium hook-and-loop closure with a leather-stiffened ring for easy adjustment by users with limited hand function.

Usable Length: 46" (117cm)  
Webbing Width: 2" (50mm)

- SH102-B2** Cinch-Mount
- SH102P-B2** Cinch-Mount (with PD306 Slip-On Pad. See pg. 8)
- SH102-B1** Flat-Mount
- SH102P-B1** Flat-Mount (with PD306 Slip-On Pad. See pg. 8)



## POWERCHAIR COMPONENTS

Improving the accessibility of power wheelchair controls lets the user achieve a higher level of independence.



## MIDLINE JOYSTICK MOUNTING SYSTEM

The Midline Joystick Mounting System improves accessibility and function. Control units and communication devices are easily mounted within reach without additional chair modifications. The precision-machined construction withstands the toughest daily use.




**MOUNTING CLAMP**  
Allows installation on the left or right armrest.  
**HW407** Tube Size:  $\varnothing$  7/8" (22mm) X  $\varnothing$  3/4" (19mm)

**TRI-LOCK™ ROTATING SHAFT**  
Operates with smooth, precision action. Prevents accidental rotation from driving position.  
**PC024**  $\varnothing$  3/4" (19mm)

**QUICK-RELEASE LEVER**  
Permits easy adjustment and removal of the Tri-Lock™ Shaft  
**FA170** For HW407 Clamp

### MIDLINE ARMS

Midline Arms support joystick boxes and communication devices weighing up to 30 lbs. (13.6 kg).

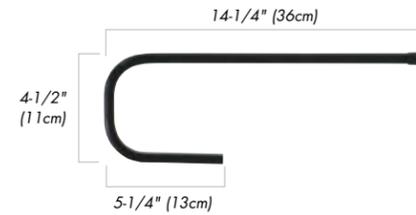

**PC031** J-Shaped Midline Arm Bracket




**PC032** L-Shaped Midline Arm Bracket




**PC033** Straight Midline Arm Bracket

### TRI-LOCK™ ROTATING SHAFT POSITIONS


Shown with a Straight Midline Arm Bracket.

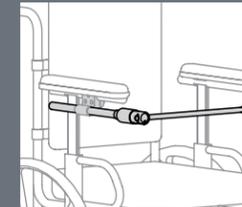


0° Horizontal driving position

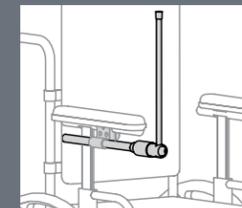

90° Allows user to transfer

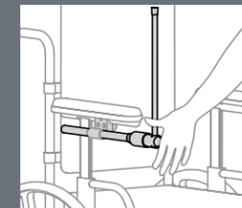

90° Press to disengage

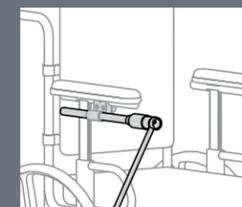

240° Moved out of the way for table access



### JOYSTICK HANDLE COLLECTION

Evaluate and select the best handle style for optimal driving control from our collection of nine flex-shaft joystick handles. Choose a display with handles to fit either 3/16" or 1/4" diameter stems.

- PC130-A KIT**  
Fits stem size  $\varnothing$  3/16"
- PC130-B KIT**  
Fits stem size  $\varnothing$  1/4"

## BODYPOINT® HAS JOYSTICK HANDLES to fit every type of user.




**3" U-SHAPED**  
**PC101A** Fits stem size:  $\varnothing$  3/16"  
**PC101B** Fits stem size:  $\varnothing$  1/4"




**4" U-SHAPED**  
**PC102A** Fits stem size:  $\varnothing$  3/16"  
**PC102B** Fits stem size:  $\varnothing$  1/4"



**BALL KNOB**  
**PC103A** Fits stem size:  $\varnothing$  3/16"  
**PC103B** Fits stem size:  $\varnothing$  1/4"

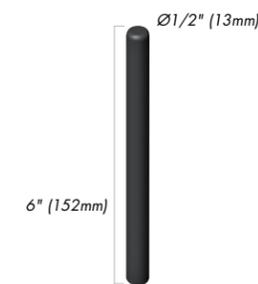

**STRAIGHT**  
**PC104A** Fits stem size:  $\varnothing$  3/16"  
**PC104B** Fits stem size:  $\varnothing$  1/4"

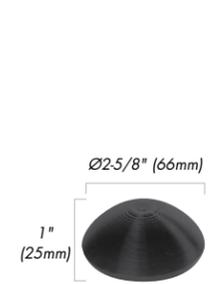

**RUBBER DOME**  
**PC107A** Fits stem size:  $\varnothing$  3/16"  
**PC107B** Fits stem size:  $\varnothing$  1/4"



**2" RUBBER MUSHROOM HEAD**  
**PC108AS** Fits stem size:  $\varnothing$  3/16"  
**PC108BS** Fits stem size:  $\varnothing$  1/4"



**2 1/2" RUBBER MUSHROOM HEAD**  
**PC108AL** Fits stem size:  $\varnothing$  3/16"  
**PC108BL** Fits stem size:  $\varnothing$  1/4"



**ALUMINUM TOP**  
**PC109A-BLUE** Fits stem size:  $\varnothing$  3/16"  
**PC109B-BLUE** Fits stem size:  $\varnothing$  1/4"  
Grey/Red available as -GREY or -RED

**MIDLINE JOYSTICK MOUNTING KITS** Kits contain: 1 Rotating Shaft, 1 Mounting Clamp, choice of Midline Arm and choice of Mounting Adapter.

KIT STYLE	JOYSTICK MOUNTING ADAPTERS				
	UNIVERSAL	INVACARE MARK IV, MPJ, PSR & PSF	QUICKIE QTRONIX	INVACARE MARK IV R2 & MARK VI SPJ*	PGDT PILOT+
J-SHAPED ARM	 PC041	 PC042	 PC043	 PC044	 PC045
L-SHAPED ARM	PC201	PC202	PC203	PC204	PC205
STRAIGHT ARM	PC201A	PC202A	PC203A	PC204A	PC205A
	PC201B	PC202B	PC203B	PC204B	PC205B

\*Needs slight modification.

## HARDWARE

Bodypoint® hardware provides a wide range of mounting options for all your seating equipment.



### FRAMESAVER™ CLAMP

Fits Tube Size: Ø7/8" (22mm) and Ø1" (25mm), Optional shim for Ø3/4" (19mm) tube also available (part no. HW303-19SHIM)

<b>HW303-B-2</b>	Quantity/Kit: 2
<b>HW303-B-10</b>	Quantity/Kit: 10
<b>HW303-B-50</b>	Quantity/Kit: 50
<b>HW303-B-100</b>	Quantity/Kit: 100



### ALUMINUM FRAME CLAMP

This double screw clamp allows you to separately secure the mounting position prior to attaching your accessory—even on an upholstered seat rail!

<b>HW310-22-A-4</b>	Ø7/8" (22mm)	Quantity: 4
<b>HW310-22-A-50</b>	Ø7/8" (22mm)	Quantity: 50
<b>HW310-25-A-4</b>	Ø1" (25mm)	Quantity: 4
<b>HW310-25-A-50</b>	Ø1" (25mm)	Quantity: 50



### QUICK-RELEASE CONNECTOR

The patented Bodypoint Quick-Release Connector provides an easily detached, solid mount that can be quickly released—even with weak hand function. The simple and reliable push-to-close, pull-to-release mechanism can be operated with one hand and locked for extra security when needed. With a slim, 1/2" profile and durable stainless steel and Delrin construction, these versatile mounts fit securely in all but the most inaccessible places. **HW311** Quantity: 1

## CLICKIT® AND FIXIT™ SEAT LATCH



Patent 5,683,201

### CLICKIT SEAT-LATCH

Spring-loaded latch release; secures solid seating drop hooks to wheelchair tubing.  
Fits Tube Size: Ø1" (25mm) and Ø7/8" (22mm) with shim.

<b>HW600-1</b>	Quantity: 1
<b>HW600-4</b>	Quantity: 4
<b>HW600-10</b>	Quantity: 10
<b>HW600-50</b>	Quantity: 50
<b>HW600-100</b>	Quantity: 100



**SEAT KIT**  
Two Clickit and two Fixit  
**HW611-1**



Patent 5,683,201

### FIXIT SEAT-LATCH

Secures solid seating drop hooks to wheelchair tubing.  
Fits Tube Size: Ø1" (25mm) and Ø7/8" (22mm) with shim.

<b>HW601-1</b>	Quantity: 1
<b>HW601-10</b>	Quantity: 10
<b>HW601-50</b>	Quantity: 50
<b>HW601-100</b>	Quantity: 100



### CAM BUCKLES Sold in pairs.

<b>FS032-2</b>	Width: 1" (25mm)
	Ø Hole: 3/16" (5mm)
<b>FS033-2</b>	Width: 1-1/2" (38mm)
	Ø Hole: 1/4" (6mm)
<b>FS034-2</b>	Width: 2" (50mm)
	Ø Hole: 1/4" (6mm)



### MOUNTING CLAMP

<b>HW404</b>	Tube size: Ø3/4" (19mm) X Ø7/8" (22mm)
<b>HW404-2</b>	Pair
<b>HW406</b>	Tube size: Ø3/4" (19mm) X 1" (25mm)
<b>HW406-2</b>	Pair

## TERMS + CONDITIONS

Bodypoint, Inc. | 558 First Ave. S., Suite 300 | Seattle, WA 98104 | [www.bodypoint.com](http://www.bodypoint.com)

### CUSTOMER SERVICE

**Toll-free: 800.547.5716**  
**Tel: 206.405.4555**  
**Fax: 206.405.4556**  
**Toll-free fax: 800.767.3828**  
**Monday-Friday, 6 AM-5 PM (PST)**

### PAYMENT INFORMATION

We accept Visa, MasterCard, or American Express

### SERVICE

Please contact Bodypoint to request a credit application, set up a new account, locate the nearest area representative or international distributor to you, or for any questions you might have. Prices, terms and conditions are subject to change without notice.

### TERMS

Standard terms are net 30 days on accounts approved for credit. 2% discount for invoices paid in full by check within 15 days of the invoice date. All shipments are FOB Seattle, Washington. Prices are subject to change without notice.

### SHIPPING

Free Domestic UPS Ground Service for orders of \$150 or more. Applies to UPS Ground Service in the U.S. Offer excludes OEM, Distributor, and VAR orders. Excludes all oversize items.

### RETURNS

Customer must contact Bodypoint within 30 days of the purchase date for a required return authorization. All returns are subject to a 15% processing fee, plus applicable freight charges.

### PAST DUE ACCOUNTS

Past due invoices will incur finance charges at the rate of 1.5% per month. Bodypoint reserves the right to hold further shipments or require prepayment for any delinquent account.

### LIMITED LIFETIME WARRANTY

All products manufactured by Bodypoint carry a limited lifetime warranty against defects in workmanship and materials arising under normal use by the original consumer. Product specifications are subject to change. For current specifications and warranty details, contact your local Independent Distributor or see our website at [www.bodypoint.com](http://www.bodypoint.com).

### LIMITATIONS ON WARRANTIES

The warranties described herein do not extend to products that have been subject to misuse, accident or tampering, nor to products that have been modified in any way other than that specifically described in the product instructions. Additional limitations may apply to distributed products not manufactured by Bodypoint—see specific manufacturer's warranty information.

Purchaser's sole remedy for defective products shall be repair, replacement or account adjustment at the discretion of Bodypoint. In all circumstances Bodypoint reserves the right to substitute an equivalent product or refund the purchase price of the original product.

The foregoing warranties are in lieu of all other warranties expressed or implied, including without limitation, implied warranties of merchantability and fitness for a particular purpose. In no event shall Bodypoint be liable for any special, indirect, incidental or consequential losses or damages resulting from the use of the product, including without limitation, loss of anticipated profits, even if Bodypoint shall have been advised of the possibility of such potential loss or damage. Bodypoint shall not be liable for any loss or damage caused by delay in furnishing product.

Bodypoint does not authorize any distributor, dealer, agent, sales person or any third party to change the terms or coverage of this warranty, and any such change shall be null and void. This warranty is non-transferable. This warranty will be governed in accordance with the laws of the State of Washington. The customer consents that jurisdiction under Venue for any action shall be in King County, Washington. This warranty information is reprinted here for information purposes only and is subject to change without notice.

### TRADEMARK INFORMATION

Aeromesh is a registered trademark of Bodypoint, Inc.  
Ankle Huggers is a registered trademark of Bodypoint, Inc.  
Bodypoint is a registered trademark of Bodypoint, Inc.  
© is a registered trademark of Bodypoint, Inc.  
Clickit is a registered trademark of Bodypoint, Inc.  
Cinch-Mount is a trademark of Bodypoint, Inc.  
Fixit is a trademark of Bodypoint, Inc.  
Flat-Mount is a trademark of Bodypoint, Inc.  
FrameSaver is a trademark of Bodypoint, Inc.  
Laminar Pad is a trademark of Bodypoint, Inc.  
Monoflex is a trademark of Bodypoint, Inc.  
Rehab Latch is a trademark of Bodypoint, Inc.  
Stayflex is a trademark of Bodypoint, Inc.  
Taperlock is a trademark of Bodypoint, Inc.

MARK IV, V and VI are trademarks of Invacare Corporation.

Qtronix is a trademark of Sunrise Medical  
Quickie is a registered trademark of Sunrise Medical  
PGDT Pilot+ is a registered trademark of PG Drives Technology, Ltd. (formerly, Penny and Giles)

© Bodypoint® 2010




558 First Avenue South, Suite 300  
Seattle, WA 98104

Presorted Standard  
US Postage  
**PAID**  
Permit # 1244  
Seattle, WA

### TRY OUR ONLINE EDUCATION

Online education allows convenient access 24 hours a day and the ability to learn at one's own pace. Our goal at Bodypoint is to provide easy to understand principles to help improve stability and function in seating with the use of primary and secondary supports. Accreditation for these courses is provided by the University of Pittsburgh-Department of Rehabilitation Science and Technology. These courses are offered free to therapists, rehab technology suppliers and other seating professionals for a limited time. Check it out at [www.bodypoint.com](http://www.bodypoint.com)

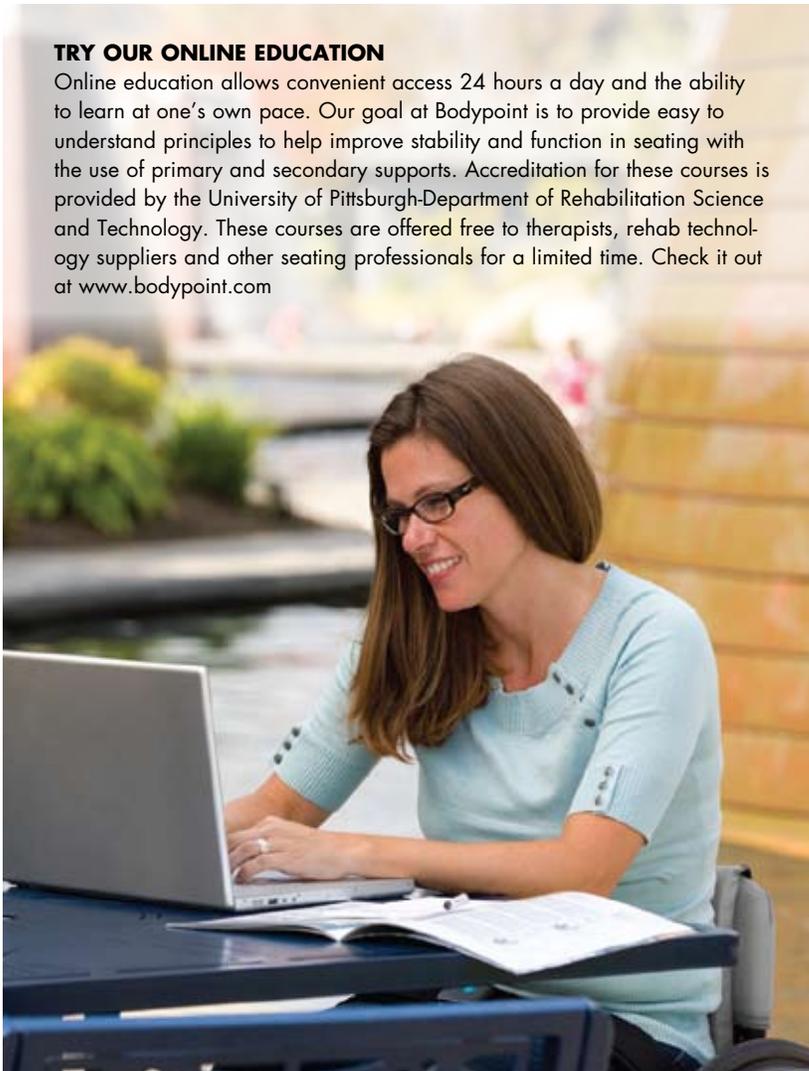

### SHOWER CHAIR BELT

The Bodypoint shower chair belt with our proprietary Aeromesh™ padding provides safety without compromising the comfort and health of the user. Water and air easily pass through the Aeromesh fabric, allowing the client to be bathed completely and the belt to dry quickly after use. It's flexible construction and universal shape allow it to be used in various ways to secure the user: across the chest or abdomen to prevent forward/lateral falls, or fitted snugly across the knees to keep clients from slipping forward and out of the chair. See page 17 for ordering information.



### PUSH-BUTTON BUCKLE

The new Bodypoint Push-Button Buckle is sleek and low-profile. Engineered for comfort and performance, this new design offers a range of features that therapists and users have requested for years. The patented construction allows for removal of the buckle cover for easy cleaning, or to change the desired button hole size—in effect providing three buckle options in a single design. See page 6 for ordering information.